



AERO TACTICAL™

TRAUMA MEDIC BERGEN



Modular &
Customisable



Optimised Visibility
& Accessibility



Built for
Endurance



Waterproof



Fits larger
CO₂ cylinder

DESIGNED WITH FORCE EXPERTISE, FOR OPERATIONAL EXCELLENCE

Elevate your state of readiness with the AeroTactical™ Trauma Medic Bergen – a meticulously crafted solution, engineered through extensive trials with various police forces across the UK. The result of this intensive product development phase is a customisable and exceptionally robust bergen, capable of accommodating modular and configurable internals. Engineered for swift frontline deployment, the **AeroTactical™ Trauma Medic Bergen** empowers confident response in every critical moment.

 **Factor of Safety 7:1**



MEDIC BERGEN: KEY BENEFITS

- ✔ Designed by forces for forces.
- ✔ Modular and configurable inserts.
- ✔ Front pouch fits any defibrillator model —→
- ✔ Hard-wearing (reinforced base and sides).
- ✔ Freestanding and upright on bottom or side.
- ✔ Zips can be security sealed.
- ✔ Chest harness increases carrying comfort.



Working Load Limit: **30 kg**


- ✔ Low light visibility with reflective strips →
- ✔ Waterproof, providing ingress protection.
- ✔ External surfaces easily wipe clean.
- ✔ Able to fit larger CO² cylinder.
- ✔ High visibility fleck on zipper pull cords.
- ✔ Document sleeve in the back panel.
- ✔ Available with 3 different coloured straps.



✔ Reinforced rubber base



✔ Reassuringly waterproof

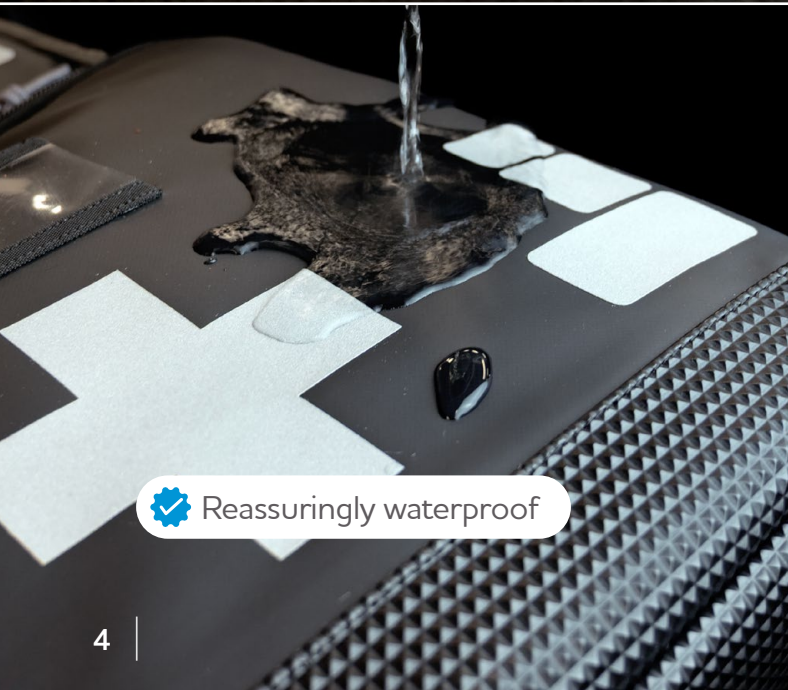

✔ VELCRO® attachment straps allows quick release of modular inserts



**EMS-PAN-BAG-BLACK****AeroTactical™ Trauma Medic Bergen – BLACK****AVAILABLE WITH DIFFERENT COLOURED STRAPS****EMS-PAN-BAG-GREEN****AeroTactical™ Trauma Medic
Bergen – GREEN****EMS-PAN-BAG-RED****AeroTactical™ Trauma Medic
Bergen – RED**

MEDIC BERGEN INSERTS: KEY BENEFITS

- ✓ Ability to insert ID tags onto pouches.
- ✓ Pouches have clear display panel to quickly identify contents.
- ✓ Quick release zips with pull cords.
- ✓ Pouch and insert surfaces easily wipe clean.



EMS-PAN-IP-L-BLUE

**AeroTactical™ Enclosed Zippered
Inner Pouch – Large BLUE**




EMS-PAN-IP-L-RED

**AeroTactical™ Enclosed Zippered
Inner Pouch – Large RED**



EMS-PAN-IP-L-BLACK

**AeroTactical™ Enclosed Zippered
Inner Pouch – Large BLACK**

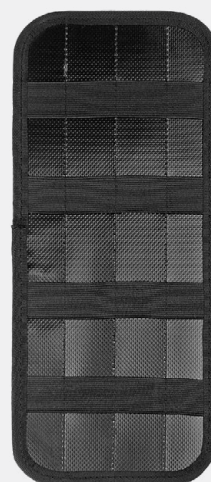

EMS-PAN-ELI

**AeroTactical™ Removeable
Elastic Loop Insert**

- ✓ High visibility lining to quickly identify any missing pouches, etc.
- ✓ VELCRO® backing allows quick release of pouch / insert.
- ✓ Waist buckle available for added comfort and security.
- ✓ First aid kit contents can be custom packed by Aero.

**EMS-PAN-IP-O2**

**AeroTactical™ O2/N2O3
Cylinder Sleeve**

**EMS-PAN-RBI**

**AeroTactical™ Removeable
Bungee Insert**

**EMS-PAN-IP-M**

**AeroTactical™ Enclosed Zippered
Inner Pouch – MEDIUM**

**EMS-PAN-IP-S**

**AeroTactical™ Enclosed Zippered
Inner Pouch – SMALL**

MEDIC BERGEN RANGE: TECHNICAL SPECIFICATIONS



AeroTactical™ Trauma Medic Bergen

Height: **59 cm**

Width: **40 cm**


Depth: **28 cm**

Internal Volume: **c. 46 L**

Working Load Limit: **30 kg**

Factor of Safety **7:1**

Proven durability: **30 kg recommended working load, SGS fatigue tested over 1,000 oscillations and separately load tested to 200 kg**

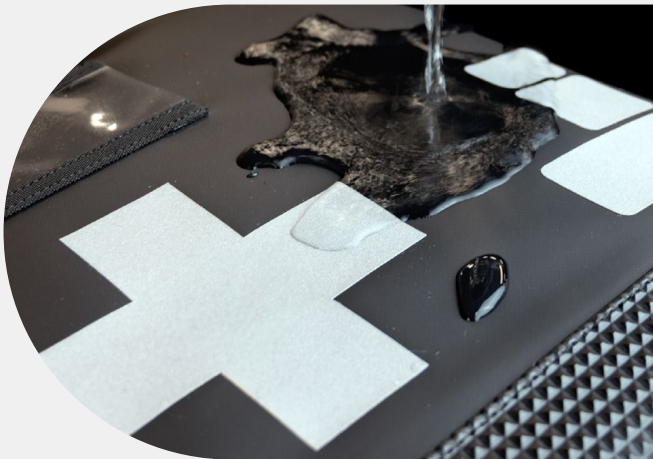

Outer Main Material:
840D Nylon with TPU coating

Interior Lining Material:
600D Polyester with PU coating

Padding (different size foams are used):
EPE foam 3mm/ 5mm/8mm; EVA foam 3mm/4mm; PE board 1mm/2mm



Webbing Material: **Polyester**

Zipper: **SBS double-runner, quick release zips with pull cords**

AeroTactical™ Enclosed Zippered Inner Pouch – Large




Height: **38 cm**

Width: **18.5 cm**

Depth: **8 cm**

Material:
500D PVC Tarpaulin

AeroTactical™ Removeable Elastic Loop Insert

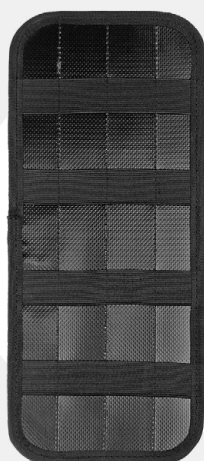

H: **33 cm**

W: **14 cm**

D: **0.5 cm**

Material:
500D PVC Tarpaulin, PE board

AeroTactical™ O2/N2O3 Cylinder Sleeve



H: **38 cm**

W/D: **12.5 cm**

Material:
PVC outer, 210D with PU coated lining, EPE foam

AeroTactical™ Removeable Bungee Insert



H: **33 cm**

W: **14 cm**

D: **1.5 cm**

Material:
500D PVC Tarpaulin, PE board

AeroTactical™ Enclosed Zippered Inner Pouch – Medium



H: **18.5 cm**

W: **18.5 cm**

D: **8 cm**

Material:
500D PVC Tarpaulin

AeroTactical™ Enclosed Zippered Inner Pouch – Small



H: **12 cm**

W: **18.5 cm**

D: **8 cm**

Material:
500D PVC Tarpaulin

Contact Aero Healthcare today

For further information on our first aid and trauma provision, defibrillators and how you could benefit, simply get in touch with a member of our team who will be happy to assist you.



HEALTHCARE | WE MAKE IT BETTER

CALL **01403 790 704**

EMAIL **sales@aerohealthcare.co.uk**

VISIT **aerohealthcare.com/uk**

© Aero Healthcare